USER MANUAL

WiFi-USB DIGITAL DEVICE

DM 93/42/CEE (MODIFIED BY DIRECTIVE 2007/47/CE) MEDICAL DEVICE IN CLASS 1

MODELS

FEDMEDMICFI

Copyright 2017, ITALECO S.R.L. All rights reserved.

Copies, whole or partial duplication, in any form, of this manual is forbidden without the written authorization of ITALECO S.R.L.

Rev. 03

Last update of this manual: Feb. 2017

WARNING! USERS SHOULD READ THIS MANUAL VERY CAREFULLY BEFORE USING THE DEVICE.

ITALECO S.R.L. won't be held responsible in case of damages or injuries due to an illegal or improper use of the product resulting from failure to follow instructions, warnings, notices or indications for use, published on this manual. "WifiMic-Fi" and "USBMic-Fi" software/app, provided together with this device, is under the Copyright. All rights are reserved. The software is protected by the copyright law and the provisions of applicable international treaties. In accordance to these laws, the licensee is permitted to use a copy of the software embedded in it in the manner provided for the operation of the products in which it is incorporated. The software or copies of the software are not for sale, all right, title and ownership of the software belongs to manufacturer or its suppliers.

For further information about the device, please contact:

MANUFACTURER

ITALECO SRL

Via Pavia 105/i, 10098, Rivoli (To), Italy

Tel: +39-011-5690297

Fax: +39-011-5690298

<http://www.mic-fi.com>

Summary

1.	UNPACKING	4
2.	PACKAGING CONTENTS	4
3.	SYMBOLS.....	5
4.	WARNINGS.....	7
5.	INTRODUCTION	8
6.	DESTINATION OF USE	8
7.	GENERAL GUIDELINES FOR HANDLING.....	8
8.	GENERAL MAINTENANCE	9
9.	INSTRUCTION FOR USE	10
9.1	DERMATOSCOPE.....	10
9.2	EARSCOPE	11
9.3	DENTALSCOPE.....	12
9.4	CAMERA for OPTICAL DEVICES	13
10.	CLEANING THE DEVICE	14
11.	DISPOSAL	15
12.	WARRANTY	16
13.	INSTALLATION AND USE OF THE SOFTWARE AND APP's	17

1. UNPACKING

We kindly remember that the packaging element (paper, cellophan, staples, tape and so on) can cut the skin or injure if handled careless. It is important to handle them carefully and by using proper equipments; do not let equipments and packaging elements near people who are not responsible.

Once you have opened the packaging, please check properly its content to be sure that all parts are inside it and that they are in perfect condition.

2. PACKAGING CONTENTS

The packaging contains:

One Mic-Fi device

One CD with software+ manual

One Charger

One USB cable

One quick guide

Accessories: depending on the intended use

3. SYMBOLS

The symbols shown on the following pages can be found on the WiFi Mic-Fi Device, on the packaging, on the box or in this manual.

WARNING SYMBOLS

WARNING This symbol shown on this manual indicates that some condition or procedures, which are not respected, could cause injuries or damages to the user or to the patient if occurs a misuse of the device.

WARNING Read the instructions

SYMBOLS FOR MAINTENANCE, TRANSPORT AND ENVIRONMENT

High

Fragile

Keep dry

CERTIFICATION SYMBOLS

Indicates that the product satisfy the essential requirements of the EU directive for medical devices 93/42/CEE

Indicates a part of B type which goes in contact with the patient.

Indicates the name of the manufacturer.

ITALECO SRL

Via Pavia 105/i,
10098, Rivoli (To),
ITALY

IDENTIFICATION SYMBOLS

Followed by an alphanumeric code indicates the lot code.

Followed by an alphanumeric code and date indicates the serial number and the production Year

OTHER SYMBOLS

Indicates that the product meet the directive 2011/65/CE related to the reduction of the use of dangerous substances in electrical and electronic devices

Indicates that the product at the end of its life must be separately collected. Thus the user must bring it to the specific collecting point dedicated to the collection of electronic and electric devices, or bring it back to the reseller in case of the purchase of a new one. The adequate separate collection contributes to a correct environment policy which allows the reuse of dismissed devices. The illegal collection will be strictly punished according to the related laws and norms.

Batteries are considered as a special waste. Do not throw them away together with domestic waste. You have to bring them to specific collection points. The illegal collection will be strictly punished according to the related laws and norms

FCC, Federal Communications Commission It is a USA commission; the symbol indicates that the product won't interfere with other devices present in your household or office

4. WARNINGS

Please read carefully these warnings before start using the Mic-Fi device. The following warnings are related to the generic use of the Mic-Fi device. A different use of Mic-Fi device could put patient's health and device's users in risk.

The device is suitable for external use only

The device does not make diagnosis

The size that can be viewed from Mic-Fi SW and App are indicative and must not to affect to the final judgment borne by the doctor.
The Mic-Fi devices are not devices with a measuring function

DO NOT INSERT THE DEVICE IN HUMAN OR ANIMAL'S ORIFICE. DANGER OF HEAVY INJURIES

Keep out of reach of children

Make sure that the covers is disinfected before and after each use

Before using the digital device please wear disposable gloves in order to keep a suitable hygienic level.

Do not direct the beam of light toward the eyes

Do not use the Mic-Fi device for capture images or videos of eyes

Do not use the Mic-Fi device on open wounds or on painful parts of the body

Do not use the device in close proximity to high heat, magnetic fields

Do not use the device in contact with lenses, glasses or other surfaces that can change the perceptions of the images

Do not use the device in WiFi mode during the charging

5. INTRODUCTION

Thank you for purchasing a Mic-Fi device. This product is designed to assist in improving the accuracy of analysis.

Only a well-documented image allows the examiners to compare their analysis with the result of other forms of investigation. In addition, a series of images derived from comparative tests repeated in a certain period of time allows an assessment of possible diseases and a more timely verification of the effectiveness of treatment prescribed or implemented to the patient. The Mic-Fi device is capable of attaching to the visual image (magnification, photos, video) a detailed verification.

Mic-Fi device despite being extremely easy to use is a medical device, and we therefore recommend reading this book before using the equipment

6. DESTINATION OF USE

The Mic-Fi device is a medical device can be observe, depends on the chosen product , capillaries, skin, and ears. The technology of this device allows you to see in real time on a common PC/tablet/smartphone. The medical device is provided with a build-in Leds lights. The medical device needs the installation of a software/application The medical device takes the Energy directly from the litium battery included and not removable

WARNING! THE DEVICE IS NOT INTENDED FOR CLINICAL INVESTIGATIONS

7. GENERAL GUIDELINES FOR HANDLING

We would always recommend wearing disposable gloves when using the device, however this can reduce the perception while operating with the device.

Do not use the Mic-Fi device with wet hands

Do not use the Mic-Fi device if damaged or if it does not work correctly

Always ensure that the device is properly connected to the hardware on which you installed the software/application supplied

Do not use the device during the charging of the battery

Do not subject the Mic-Fi device to sudden shocks

8. GENERAL MAINTENANCE

Check the product and accessories before use, even after it has been cleaned. Disinfection of the instrument is needed before use (the assessment has to be performed by the doctor). The product is made of materials resistant to corrosion and environmental conditions provided for normal use, and does not require specific attention, but it must be kept in a closed environment, clean and dry, avoiding exposure to light and protecting it from dust in order to ensure the hygiene conditions.

Do not place the device and accessories in direct sunlight, high temperatures or in humid or dusty places

After using the device always separate from accessories and specula, cleaning them with a mild disinfectant solution, then proceed with the procedures for cleaning and disinfection the rest of the equipment, put headsets and specula in the case provided, taking care to disconnect the Mic-Fi from power supply or USB port

Do not immerse the device in liquids

DO NOT EVER DISASSEMBLE THE DEVICE

9. INSTRUCTION FOR USE

9.1 DERMATOSCOPE

- a) Put the part to be observed below and in the centre of the illuminated area of your Mic-Fi device
- b) Turn the wheel of focus to maximise viewing and have a sharper image
- c) Click on "capture" from software/App to take the picture that interests you

Note: For optimal use of the device without support, grip the device between thumb and index accommodating it in the palm of the hand and using your thumb to rotate the wheel of focus. The tool can be rotated on its axis in order to achieve the rotation of the image on the screen.

Please ask to the patient to assume a comfortable position and not to move during the whole visit.

- Connect the device to the computer/tablet or smartphone
- Launch the Mic-Fi software (previously installed) See paragraph 13
- Place the device on the part to be observed
- Turn the dial to get a clear picture
- Choose the "capture" button to get a photo or video.
- Move the device from the patient
- Tell the patient is able to move again

The device is intended for external use only Do not enter the device in the oral cavity

Do not use the device in contact with wounds or painful parts of the body.

The size that can be viewed from Mic-Fi SW and App are indicative and must not to affect to the final judgment borne by the doctor.

9.2 EARSCOPE

This medical device is used in conjunction with a speculum (sizes 3, 4, 5, 6, 7mm) which allows you to perform visual examinations of the hearing canal outside of the walls of the eardrum through direct analysis of otoscopic photos or digital video caught by the device. In this configuration the device offers doctors a valuable method for improving the accuracy of the visual analysis.

Mic-Fi otoscope requires the deployment of a speculum CE marked disposable. Put the speculum on the device before capturing pictures or video (simple joint) and ask the patient to take a comfortable position and to remain immobile for the duration of the visit

- Connect the device to the computer/tablet or smartphone
- Launch the Mic-Fi software (previously installed) See paragraph 13
- Place the device on the part to be observed
- Turn the dial to get a clear picture
- Choose the “capture” button to get a photo or video
- Move the device from the patient
- Tell the patient is able to move again

Earscopic investigations must be achieved by keeping the device parallel to the floor of the room in order to minimize the risk of eardrum perforation caused by the fall of the device or speculum

The device is intended for external use only. Visual surveys of the ear and the hearing canal must be conducted exclusively by using a speculum (CE marked)

The size that can be viewed from Mic-Fi SW and App are indicative and must not to affect to the final judgment borne by the doctor.

9.3 DENTALSCOPE

In the dental field the device with the appropriate accessories holder for dental mirror may be used to perform visual examinations of the oral cavity, through direct analysis in real time on the monitor or later through the analysis of the video and images captured and stored via software. In this configuration the device facilitates the dentist during the inspection of the oral cavity allowing a better view and a better working posture and ergonomics. Put the mirror into the holder device and ask the patient to take a comfortable position and to remain immobile for the duration of the visit

- Connect the device to the computer/tablet or smartphone
- Launch the Mic-Fi software (previously installed) See paragraph 13
- Place the device on the part to be observed
- Turn the dial to get a clear picture
- Choose the “capture” button to get a photo or video
- Move the device from the patient
- Tell the patient is able to move again

The device is intended for external use only Do not enter the device in the oral cavity

Make sure that the holder has been regularly cleaned and disinfected before use

Do not to use the device simultaneously with the cleaning and dental descaling

The size that can be viewed from Mic-Fi SW and App are indicative and must not to affect to the final judgment borne by the doctor.

9.4 CAMERA for OPTICAL DEVICES

Connect the Mic-Fi medical device to your optical instrument
The connection can be direct (through the step C) or by C-coupler
Visit www.mic-fi.com for viewing various accessories

- Connect the device to the computer/tablet or smartphone
- Launch the Mic-Fi software (previously installed) See paragraph 13
- Place the device on the part to be observed
- Choose the “capture” button to get a photo or video

The size that can be viewed from Mic-Fi SW and App are indicative and must not to affect to the final judgment borne by the doctor.

10. CLEANING THE DEVICE

It is very important to avoid touching the lens of the device because it is very delicate. The lenses should be cleaned regularly, as particles of dust on them could have a negative impact on image quality. For this reason, it is recommended to protect the device with its protective casing when not in use. Avoid using any other objects or liquids that could easily scratch or damage the surface of the lens or the sensor. To clean the body of the device use a slightly dampened soft cloth. The product is designed to last a lifetime, and is resistant to most chemicals products except ketones MEK, solvents and substances used to remove the smalt.

The cleaning and disinfection operations should take place after each use.
DO NOT USE the device for at least 15 minutes after its cleaning operation

Avoid touching the surface of the lenses When cleaning the lens, if necessary, use dry compressed air to remove dirt

Do not use abrasive materials to clean the lens, it could cause permanent damage to it

11. DISPOSAL

The crossed waste bin indicates that the product at the end of its useful life must be collected separately from other waste.

Indicates that the product at the end of its life must be separately collected. Thus the user must bring it to the specific collecting point dedicated to the collection of electronic and electric devices, or bring it back to the reseller in case of the purchase of a new one. The adequate separate collection contributes to a correct environment policy which allows the reuse of dismissed devices. The illegal collection will be strictly punished according to the related laws and norms.

HOW TO TREAT BATTERIES AND ENERGY CHARGER.

According to the articles 9 and 22 of the Italian Dlgs of the 20th November 2008, n°188 "Implementation of the directive 2006/66/CE related to batteries and battery accumulators which repeal the directive 91/157/CE". The collection and recycling of batteries contributes in reducing the potential environmental impact on health created by the substances held in these components
ATTENTION! Batteries are strictly considered as special wastes! Do not put them together with normal wastes. Users have to collect them in a special way or bring them to special collection points. Every shop in which batteries are sold have to dispose of one collection point. And the collection of batteries is not subordinate to the purchase of new batteries. The symbol placed on batteries and accumulators indicates that batteries must be separately disposed from household wastes, according to EU Directive 2006/66/CE and with the local laws. If under the symbol there is also a chemical symbol, in oservance to the directive, this symbol indicates the presence of metal (Hg = Mercury, Cd = Cadmium, Pb = Lead) in the battery with a higher level than the one allowed by the directive. The illegal disposal of batteries will entail the application of administrative sanction, according to D. Lgs. n. 22/1997" (article 50 and following of D.Lgs.n.22/1997)

12. WARRANTY

Congratulations for purchasing our product. This product HAS BEEN MANUFACTURED TO THE HIGHEST QUALITY AND STANDARDS. Italeco guarantees this product against any defects due to materials or workmanship for a period of 2 years from the date of original purchase. To take advantage of the guarantee, you must retain your purchase invoice. If the product is found to be defective within the warranty period it should be returned to the selling agent with proof of purchase. The warranty is null and void if the damage is resulting from carelessness or misuse. No compensation will be given by Italeco and its distributors if the product is returned during the warranty period. In addition, any replacement or repair carried out during the warranty period, will not prolong the duration of the guarantee.

The warranty lapses in the following cases:

- The product has been tampered with or the label placed on the usb cable bearing the serial number has been removed, damaged or rendered unreadable
- The product has been damaged or has been misused
- The product has been opened or repaired by unauthorised personnel
- The product has been damaged due to non compliance with the instructions in the booklet

Products considered to be defective should be returned only to the dealer from which it was purchased accompanied by the invoice.

The product will be repaired or replaced, under our sole discretion.

ITALECO SRL

Via Pavia 105/i, 10098, Rivoli (To), Italy


Tel: +39-011-5690297

Fax: +39-011-5690298

13. INSTALLATION AND USE OF THE SOFTWARE AND APP's

13	INSTALLATION AND USE OF THE SOFTWARE AND APP's	17
13.1	Use of the device	18
13.2	Charging Instructions.....	18
13.3	Working at Wi-Fi Mode	19
13.4	Setting for Internet and Microscope's Wi-Fi Connection at the same time.....	19
13.5	Wi-Fi Security Setup	20
13.6	"Mic-Fi" on iPad/iPhone	21
13.6.1	Installation on iPhone / iPad	21
13.6.2	Calibration on iPhone / iPad	22
13.7	"Mic-Fi" on Android Tablet /Phone	23
13.7.1	Installation on Android Tablet / SmartPhone.....	23
13.7.2	Calibration on Android Table / Phone	24
13.8	"WiFiMic-Fi" on Windows Computer	25
13.8.1	File	26
13.8.2	Setting.....	26
13.8.3	Video Processing.....	28
13.8.4	Images on Windows	28
13.8.5	Calibration on Windows Systems	29
13.9	"Mic-Fi" on Mac Computer.....	31
13.9.1	Installation on Mac Computer	31
13.9.2	Images on MAC.....	32
13.9.3	Calibration on Mac Computer	33
13.10	Working at Wired USBMode	36
13.11	"USBMic-Fi" on Windows	36
13.11.1	Running "USBMic-Fi"	36
13.11.2	Menu.....	37
13.11.3	Setting.....	38
13.11.4	Calibration on Windows	39
13.12	"USBMic-Fi" on MAC OS	40
13.12.1	Running "USBMic-Fi"	40
13.12.2	Setting.....	41
13.12.3	Calibration	41


10.1 Use of the device

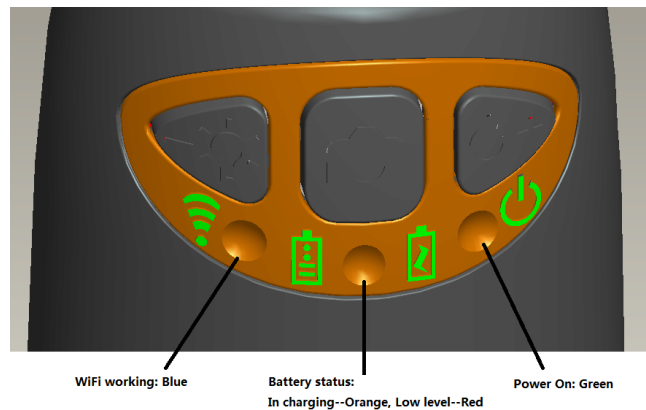- **Wi-Fi:** Press  for 2 seconds; the microscope will generate a Wi-Fi network
- **Wi-Fi Indicator**  (**Blue LED**): Flashing light indicates data transmission; Constant light indicates no data transmission.
- **Power Indicator**  (**Green LED**): Lights when power on Wi-Fi or connected to PC or AC adaptor via MiniUSB cable.
- **Led :** Press  for 2 seconds to switch the led; press + and – to increase or decrease the intensity of light.
P.S. this key does not work in the FEDMEDMICFIEYE/EYE5M models
- **Photo:** Press  to capture the picture directly from microscope
- **Snapshot Button:** Every time after the success connection of the WiFi microscope, you have to enable the Snapshot feature on App (The LED illumination will be off and then lights up again). The Snapshot only works with the receiver device which the snapshot feature of its app is enabled in the last.

10.2 Charging Instructions

1. Connect the DC Adapter (plug the AC adapter into an AC outlet) or PC's USB cable, it charges whether the power switch is on Wi-Fi or OFF/USB. **We DO NOT recommend charging Wi-Fi Microscope while using**
2. **Battery led**  (**Orange LED** or **Red LED**): **Red LED** when the battery level is low; **Orange LED** when the battery is charging
3. **The charge led** : When the battery is at full charge, the battery LED comes OFF

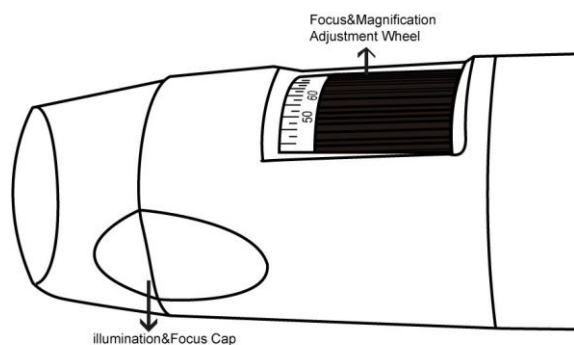WARNING : LACK OF IMAGE OR SLOW VIDEO ARE THE RESULT OF A LOW BATTERY. BEFORE THE FIRST USE OF THE MICROSCOPE, PLEASE MAKE A CHARGE OF AT LEAST 12 HOURS.

10.3 Working at Wi-Fi Mode

Before start Mic-Fi app, please make sure the Wi-Fi microscope has already connect to your display device by Wi-Fi.

1. Install:
 - a) For **iPad** and **iPhone**: Search “**Mic-Fi**” in **App Store** to download and install it.
 - b) For **Smart-Phone** and **Tablet Android**: Search “**Mic-Fi**” in **Google Play** to download and install it.
 - c) For **Windows** and **Mac PC**: Download the Mic-Fi Software by www.mic-fi.com
2. Press for 2 seconds , wait for ~20 seconds until the Wi-Fi indicator  and the illumination LEDs light, then the Wi-Fi signal will be transmitting correctly.
3. Check the wireless connection on your Phone/Tablet or PC and a list of available networks will appear. Select the SSID “Wi-Fi2SCOPE_xxxxxx” network from the list to connect; the network “Wi-Fi2SCOPE_xxxxxx” requires a password: enter the password “**12345678**” and click Ok.
4. Once connected, please wait for about 5 seconds to open the “Mic-Fi” application to see your images. If there is no video display, please close and restart “Mic-Fi” app again. The Wi-Fi indicator LED  flashes and video stream signal should be displayed in the app. Video streaming can be broadcasted to many devices, but the frame rate of live video will decrease as the number of receiving devices increase.


10.4 Setting for Internet and Microscope's Wi-Fi Connection at the same time

You can use the Wi-Fi network of the microscope and that of the PC / SmartPhone / Tablet at the same time

1. To access the internet by Wi-Fi connection at the same time the Wi-Fi microscope is working, open an internet browser, input <http://10.10.1.1> and click “enter”. A pop-up window will require you to enter a username and password, please enter: **admin** and **admin**, and then this window should appear:


2. Click  in the window to pop-up the screen:

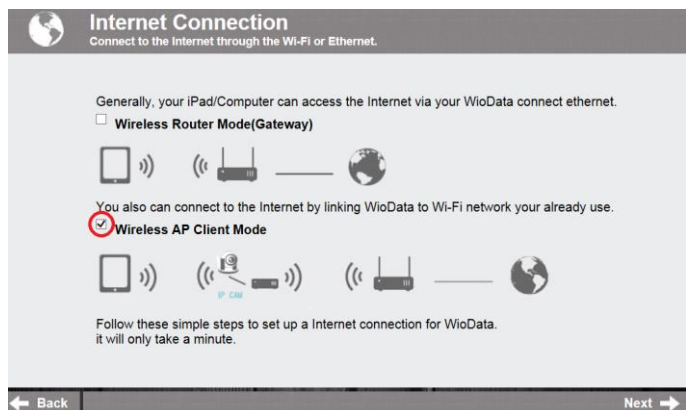3. Select “Wireless AP Client Mode”, click **Next** to pop-up the screen:

4. Click **Refresh** till your desired network SSID appears, and click it to pop-up the screen:

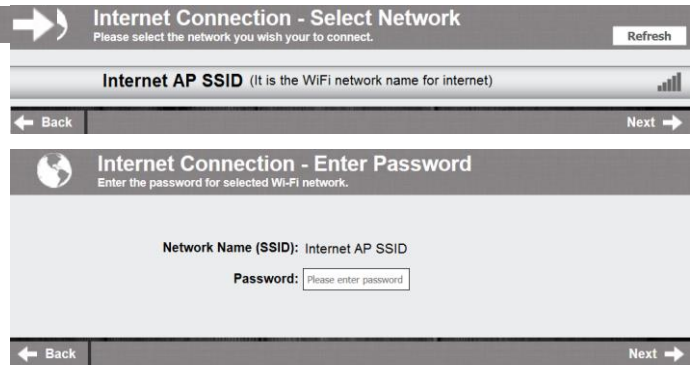Internet Connection - Select Network
Please select the network you wish to connect.

Internet AP SSID (It is the WiFi network name for internet)

← Back Next →

Internet Connection - Enter Password
Enter the password for selected Wi-Fi network.

Network Name (SSID): Internet AP SSID

Password:

← Back Next →


10.5 Wi-Fi Security Setup

To change the Wi-Fi access password:

1. Open an internet browser, input <http://10.10.1.1> and click “enter”. A pop-up window will require you to enter a username and password, please enter: **admin** and **admin**, and then this window should appear:

2. Click  in the window to pop-up the screen:

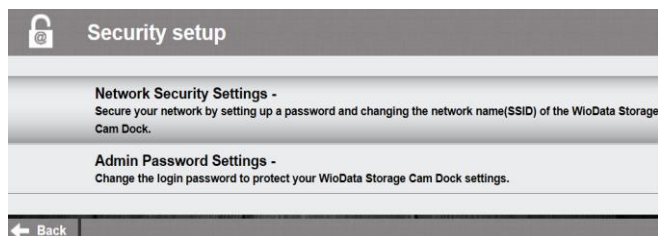
Security setup

Network Security Settings -
Secure your network by setting up a password and changing the network name(SSID) of the WioData Storage Cam Dock.

Admin Password Settings -
Change the login password to protect your WioData Storage Cam Dock settings.

← Back

3. Click “Network Security Settings” to pop-up the screen on the right.
4. Enter your desired Network Name (SSID) and password in “WPA-PSK Key” and “Re-Enter”, and click “Next”. The setting is done, now you can re-connect the new SSID Wi-Fi by entering the new password.

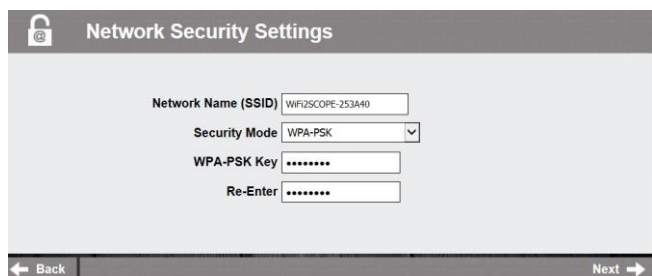Network Security Settings

Network Name (SSID)

Security Mode

WPA-PSK Key

Re-Enter

← Back Next →


10.6 "Mic-Fi" on iPad/iPhone

10.6.1 Installation on iPhone / iPad

Search "Mic-Fi" in App Store to download and install it.

Click the  icon to display video.

Warm-up Tips: If the live video is black or slow, please restart "Mic-Fi" app again.

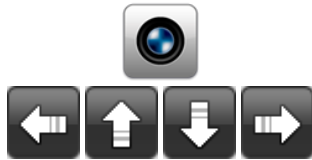Capture photo or record video (The pixel resolution has to be set at 640x480 when recording video).

Navigation icon - moving a dot left/up/down/right for calibration function

Photos folder or videos folder.

Calibration.


Switch between photo taking and video recording.

Resolution Setting. (When changing the resolution, please wait until the illumination LED is off and lights up again – it takes about 5 seconds).

Resolution setting and Default:

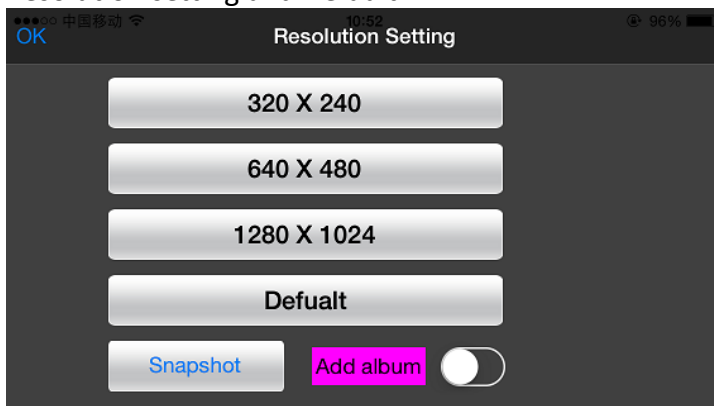: Enable the button snapshot.


: Add new album folder (Default is OFF)

10.6.2 Calibration on iPhone / iPad

Calibrate the device at every magnification


The size that can be viewed from Mic-Fi SW and App are indicative and must not to affect to the final judgment borne by the doctor.

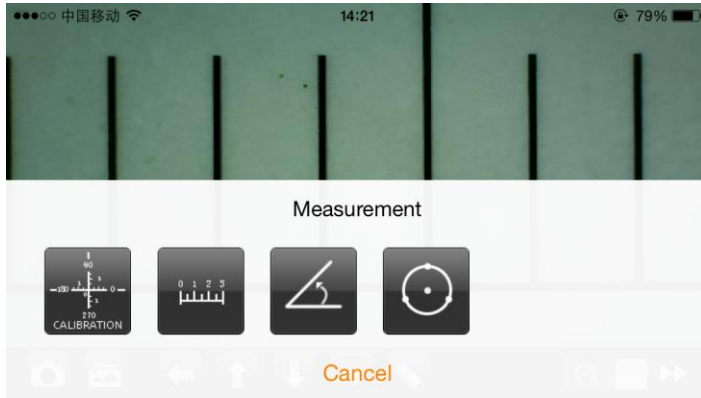
Use the enclosed calibration ruler as to-be-observed object. Bring the image into focus by turning the "Focus & Magnification adjustment wheel";

Click  to pop-up the screen on the left.

Click "Calibration"  to back the live video, click two dots (moving the dot by clicking  arrows which draw a your desired line, then click on the live video to show the following screen:

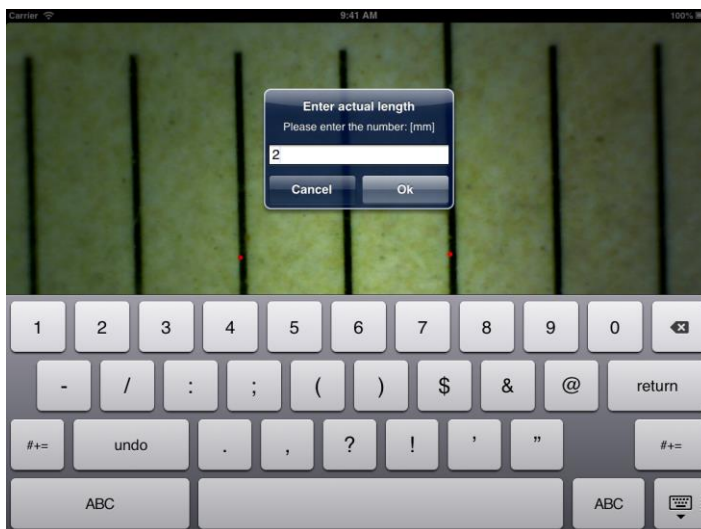Enter the value of length and press "OK"

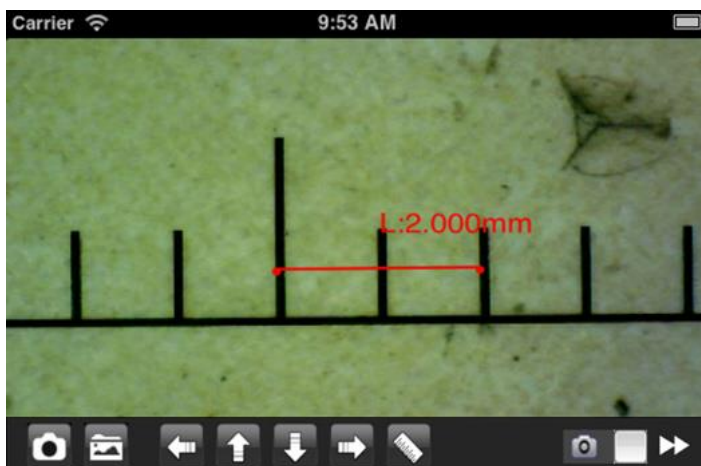Now the calibration is completed. (Don't adjust again the "Focus & Magnification wheel"). If the new photo object is not in focus, please only adjust the distance between the microscope and the object until the photo is clear.

Line Allows you to view the linear distance between two points.
Simply click and drag to the desired length, and click again to finish.

Angle Allows the viewing of an angle. Start at the pivot point and extend out to start measuring an angle.


Circle Allows the viewing of the circumference, area, and radius of a circle.

Click the ruler icon  and select “line”.

Select on the image the two points whose distance is to be viewed and press again on the screen.


10.7 “Mic-Fi” on Android Tablet /Phone

10.7.1 Installation on Android Tablet / SmartPhone

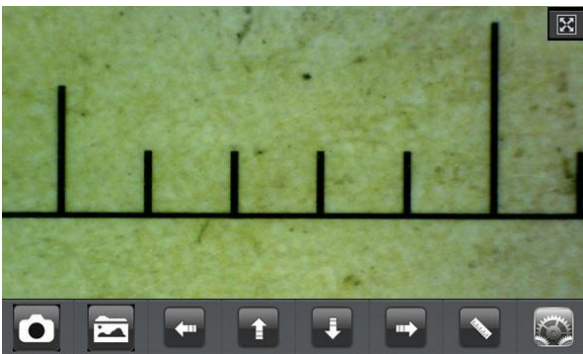Install the “Mic-Fi” app and run it by click the icon 

If the live video is black or slow, please restart “Mic-Fi” app again.

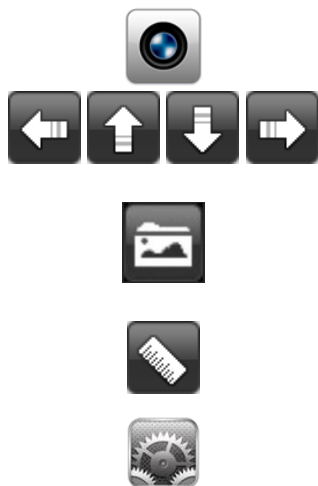Capture photo.

Navigation icon—moving a dot left/up/down/right for calibration operation.

Photos folder.

Calibration

Resolution Setting. (When changing the resolution, please wait until the illumination LED is off and lights up again – it takes about 5 seconds).

Resolution Setting(3.7)

320x240	☐
640x480	☒
1280x1024	☐
Snapshot	Rename(V)
Default	

Snapshot

: Enable the Button Snapshot.

Rename(V)

: Enable to rename the photo for every capturing

Default


: Default.

10.7.2 Calibration on Android Table / Phone

Calibrate the device at every magnification

The size that can be viewed from Mic-Fi SW and App are indicative and must not to affect to the final judgment borne by the doctor.

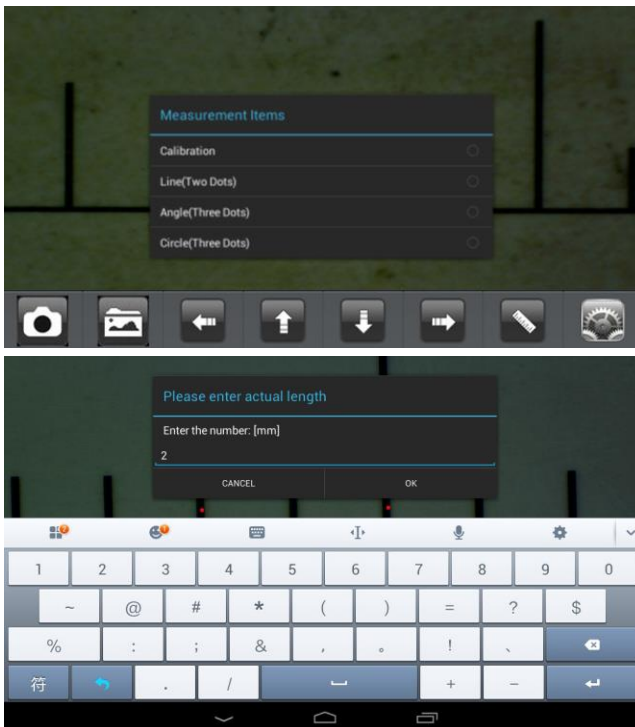Use the enclosed calibration ruler as to-be-observed object.


Bring photo into focus by turning the “Focus & Magnification wheel”; Click on  to show the screen on the left.

Click “Calibration” to come back to the live video.

1. Click to select a point and refine the selection by clicking the arrows

2. select in the same way a second point at a known distance; in this way will be shown with a line of known length
3. press any point on the screen
4. Enter the value of length and press "ok"

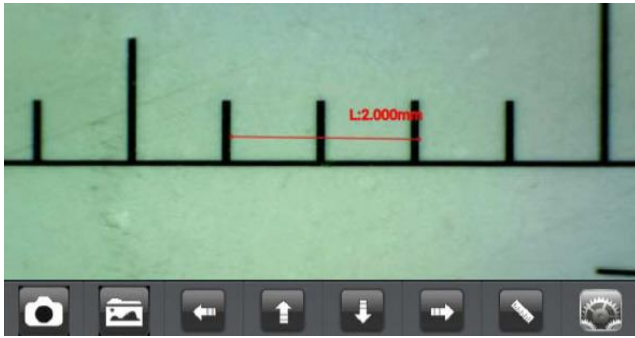Now the calibration is completed. Click on screen again (**Please, do not adjust again the “Focus & Magnification wheel”**). If the new photo object is not in focus, please only adjust the distance between the microscope and the object until the photo is clear.


Line Allows you to view the linear distance between two points.
Simply click and drag to the desired length, and click again to finish.

Angle Allows the viewing of an angle.

Circle Allows the viewing of the circumference, area, and radius of a circle.
Click the ruler icon  and select “Circle”.
Select on the image three points on the circle and press again on the screen.


Click  to save the image (**Please, do not adjust again the “Focus & Magnification wheel”**). If the new photo object is not in focus, please only adjust the distance between the microscope and the object until the photo is clear.

13.8 “WiFiMic-Fi” on Windows Computer

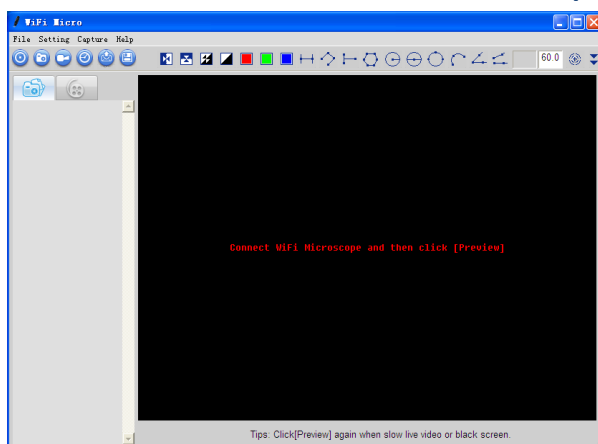
1. After the installation of “WiFiMic-Fi” the icon  will appear. Click on the icon to run it. The download is available from www.mic-fi.com

Click “Preview” icon and the Mic-Fi main screen will display video as below: (**it waits about 5 seconds to display video; if the live video is black or slow, please click “Preview” icon once again to try**).

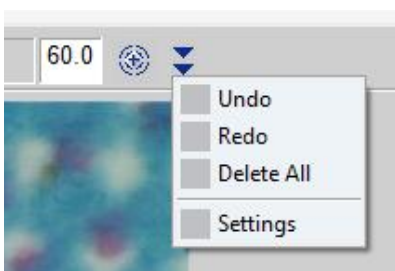-  Preview, live image
-  Take a picture
-  Take/Stop Video
-  Take time-lapsed pictures
-  Email the selected photos
-  Save the selected photo or video file

13.8.1 File

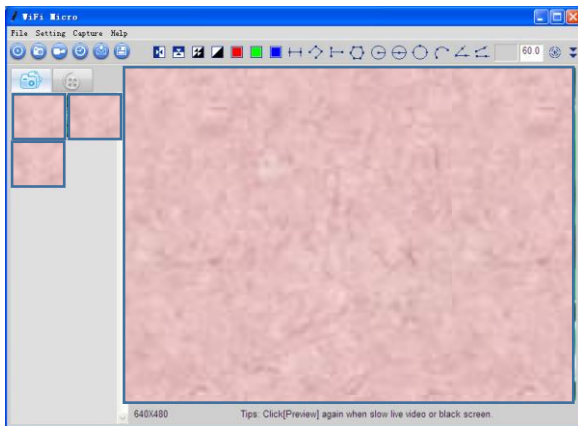
Select a photo or video and click “File” on the menu bar.

- **Open:** Open the selected photo or video to be displayed 
- **Save as** : Save the selected photo or video file as the desired folder and file format.
- **Delete:** Delete the selected photo or video.
- **Delete all:** Delete the all photos or videos in the photo or video folder.
- **E-Mail**  : Email the selected photos using your default email account.
- **Favorites:** Open the folder containing the captured photos and recorded videos.


13.8.2 Setting

Video Setting

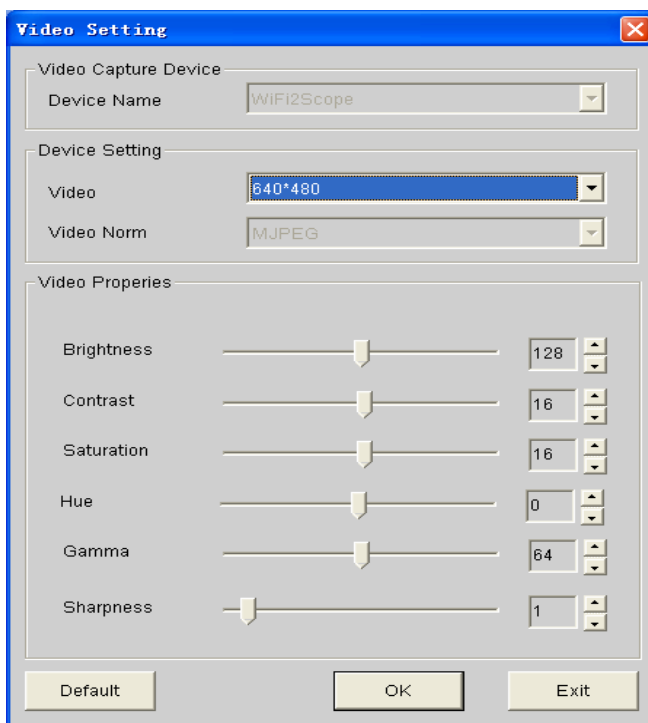When you change the resolution, please wait until the illumination LED is off and lights up again before clicking OK - it takes about 5seconds.

Setting

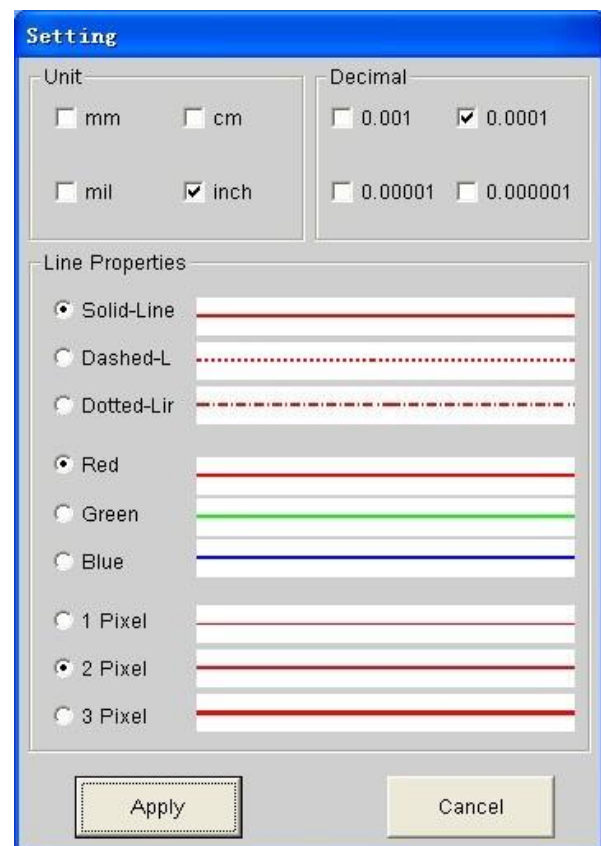Change language: The language selection allows you to select the WiFiMic-Fi interface language.

Snapshot Enable: Enable the button snapshot. The video will be stop, please Click “Preview”  again.

Warm Tips: If the Windows Security Alert appears when enable snapshot (Windows 8 as example) as below window:


You have to select the both check boxes of “Private networks” and “Public networks”, and then click “Allow access” to allow “Snapshot” program communicate on these networks, or the button snapshot doesn’t work. If miss the Allow operation, please click “Control panel->Windows firewall->Restore defaults”, and restart Mic-Fi and “Snapshot Enable” to reset it.

Language: The language selection allows you to select the Mic-Fi interface language.

Take Picture:  Capture a photo


Take/Stop Video  In the active live video window, click to start/stop recording

Recorder Setting

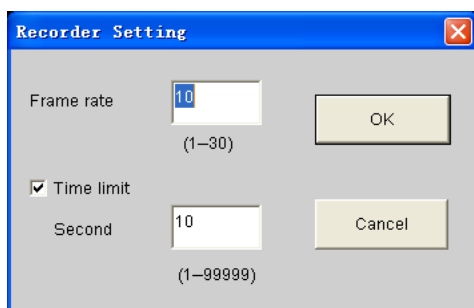- **Frame rate:** is the number of frames per second that are displayed. The number range is 1~30
- **Timer limit:** is the time of video recording to avoid video file size exceeding the disk space. The number range is 1~99999
- **Take time-elapsd pictures:** Press the time-lapsed button  on the active preview window

Capture Setting

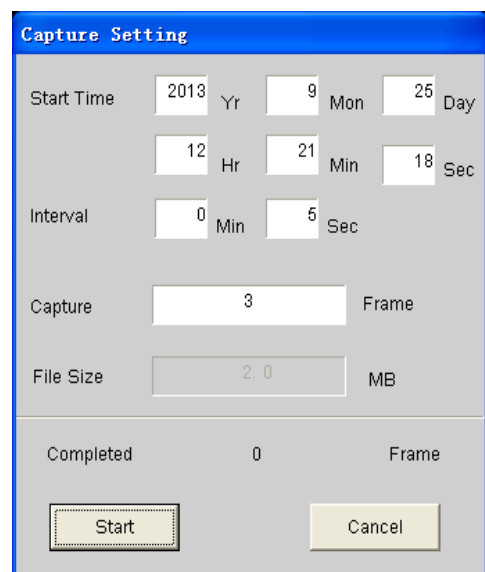13.8.3 Video Processing

Please click “Preview”  to resume the live video after finished operation of video processing below.


-  Horizontal Mirror
-  Vertical Mirror
-  Negative
-  Binary Segmentation
-  Monochrome display (Red)
-  Monochrome display (Green)
-  Monochrome display (Blue)
-  Calibration on Video

13.8.4 Images on Windows

You can double-click the selected Photo and view it as below shown.

Photo file name: yyyyymmddhhmmss_nnnx_Nxxxx is a combination of date, time, magnification and serial number. For example,

Still_2009-03-06_102949_060X_N0002 figures: Number is 0002, magnification is 60, date is Mar-06-2009, time is AM10:29:49.

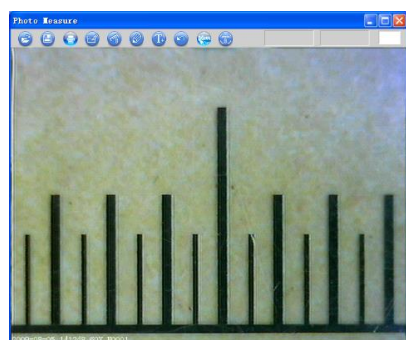Open a photo
from a folder

Save as

Print this photo

Select a target

Measure function items

Pen setting

Add text into photo

Delete the selected
target

Move the selected
target


Calibration

13.8.5 Calibration on Windows Systems


Calibrate the device at every magnification


The size that can be viewed from Mic-Fi SW and App are indicative and must not to affect to the final judgment borne by the doctor.

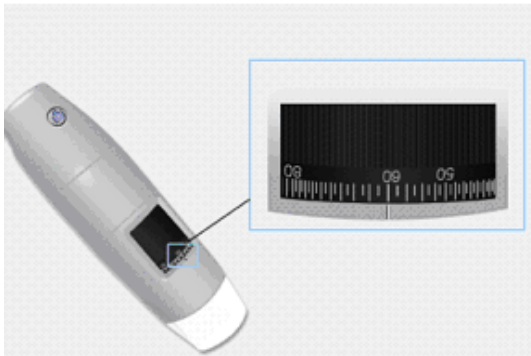
1. To calibrate, select a known measurement standard, such as a ruler, as to-be-observed object.
2. Bring the image of the selected ruler into focus by turning the "Focus & Magnification wheel"
3. Determine the magnification by taking a reading from the scale on the focus dial as below (The default magnification is 60x)

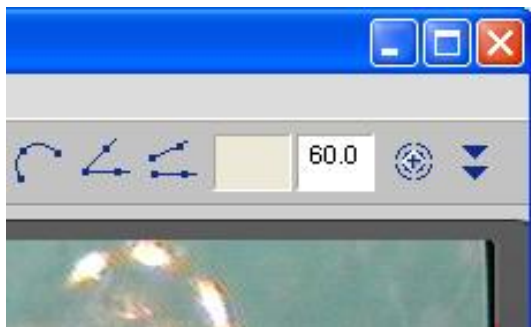4. Enter the scale value in the "Scale" box in the right-hand scale box of the screen

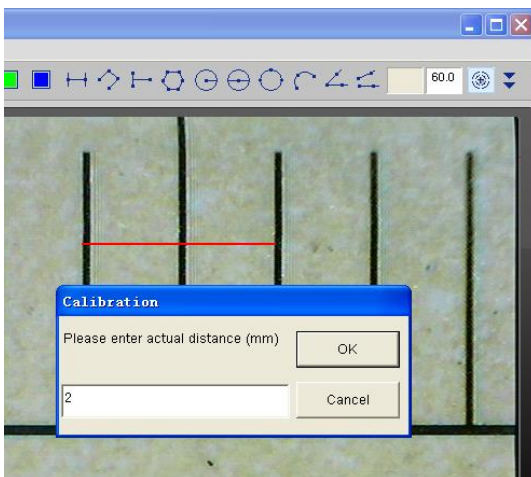5. In the live video window click on the calibration button , then click and drag to 2mm length, a small window pop up
6. Enter the actual length (the example is 2mm) and click "OK" button.

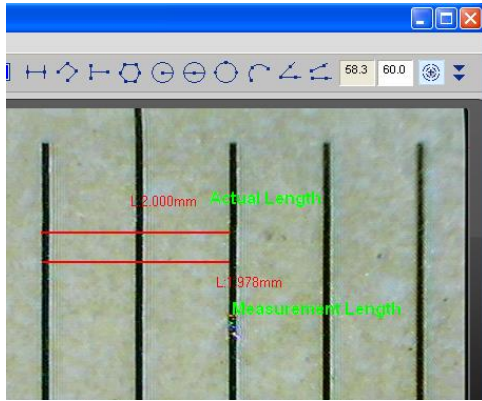The software will automatically calibrate the scale of magnification and will show the calibrated scale of magnification at the right-hand box of calibration button

Now the calibration is completed (Please, do not adjust again the "Focus & Magnification wheel").

If the new photo object is not in focus, please only adjust the distance between the microscope and the object until the photo is clear.

Line

Allows you to view the linear distance between two points. Simply click and drag to the desired length, and click again to finish.

Point to Line

Allows you to view of a line 90 degree from a certain line. Simply create a line to represent the base by clicking once to start, drag, and then click again to set the endpoint.

Continuous Line

Allows you to view of the distance between multiple connecting lines. Simply click and drag to form one section of distance, click again to start another section. Continue until the total desired distance is viewed. Double-click to finish.

Polygon

Allows you to view of a polygon. Simply click and drag to form the desired length, and click again to start the next section, double-click to finish the polygon viewed.

Radius of Circle

Allows you to view the circumference, area, and radius of a circle. Simply click and extend out to the desired radius.

Diameter of Circle

Allows you to view the circumference, area, and radius of a circle. Simply click and extend out to the desired diameter.

Three Point Circle

Allows to view of the circumference, area, and radius of a circle. Simply click on any three points on the circle you wish to view.

Three Point Arc

Allows to view of an arc. Click on three consecutive points on an arc to view.

Three Point Angle

Allows to view of an angle. Start at the pivot point and extend out to start viewing an angle.

Four Points Angle

Allows to view of an angle by selecting four points. Simply select two points from one line segment and another two points from another line segment to complete the angle viewing.

13.9 “Mic-Fi” on Mac Computer

13.9.1 Installation on Mac Computer

Installing the “Mic-Fi” in the OS X 10.8 or later version, the system will pop-up the message as below window

You need to set your computer as below:


Click Apple icon  to pop-up the window on the right, click “System Preference” to show the following window.


Click “Security&Privacy” to show the following window.

Click “Click the lock to make changes” to pop-up the following window.

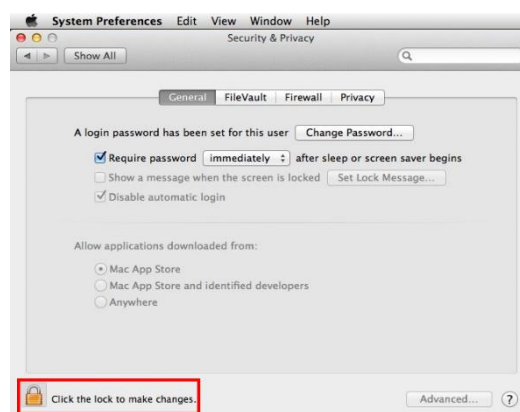Input “Name” and “password” of your system, and click “Unlock” to show the following window:

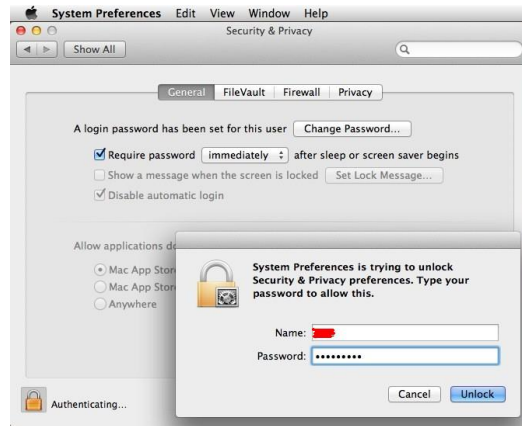Select “Anywhere” to finish the setting.
Now you can install “Mic-Fi” Software.


Double-click “WiFi Mic-fi” icon  after completing installation.

Click “Preview”  and then the WiFiMic-Fi main screen will display video (it waits about 5 seconds to display the live video; if the video is blank or slow, please click “Preview” icon once again).

13.9.2 Images on MAC

-  Preview (Connect the Wi-Fi device and display video)
-  Capture photo
-  Photo folder: Open the folder which contains the captured photos
-  Setting (Show / Hide)
-  “About”

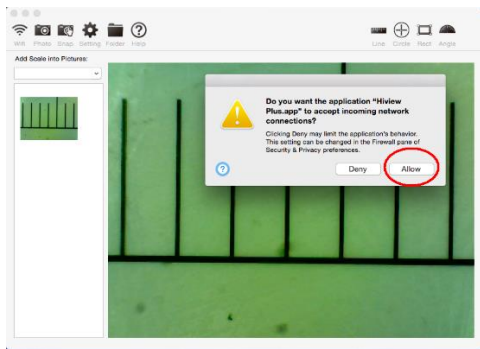Enable the Button Snapshot feature.

When change any item of the Setting or enable the button snapshot, you have to wait until the illumination LED is off and lights up again - it takes about 5seconds.


13.9.3 Calibration on Mac Computer

Calibrate the device at every magnification


The size that can be viewed from Mic-Fi SW and App are indicative and must not to affect to the final judgment borne by the doctor.


1. Select a known measurement standard as the target object, such as a ruler, to be used for calibration. Bring the image of the selected ruler into focus by turning the “Focus & Magnification adjustment wheel”, and determine the magnification by taking a reading from the scale on the focus dial. Click  icon button to capture an image and then double-click the image, One window appear as below:

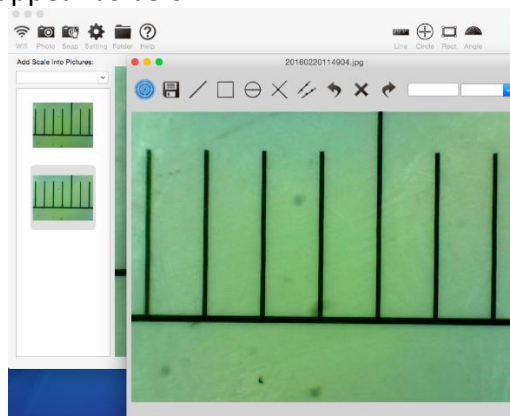
2. Click  button, the calibration window appear and enter; Title: the calibration name, such as “Calibration1”; Mag: the magnification which read from the microscope wheel, such as “60”; Length: the actual length which you will draw line, such as “2” . Finally, Click “Add Calibration” when you finish above operation, and then directly draw a line **(do not click “line” button)** which length is 2mm to complete the calibration.

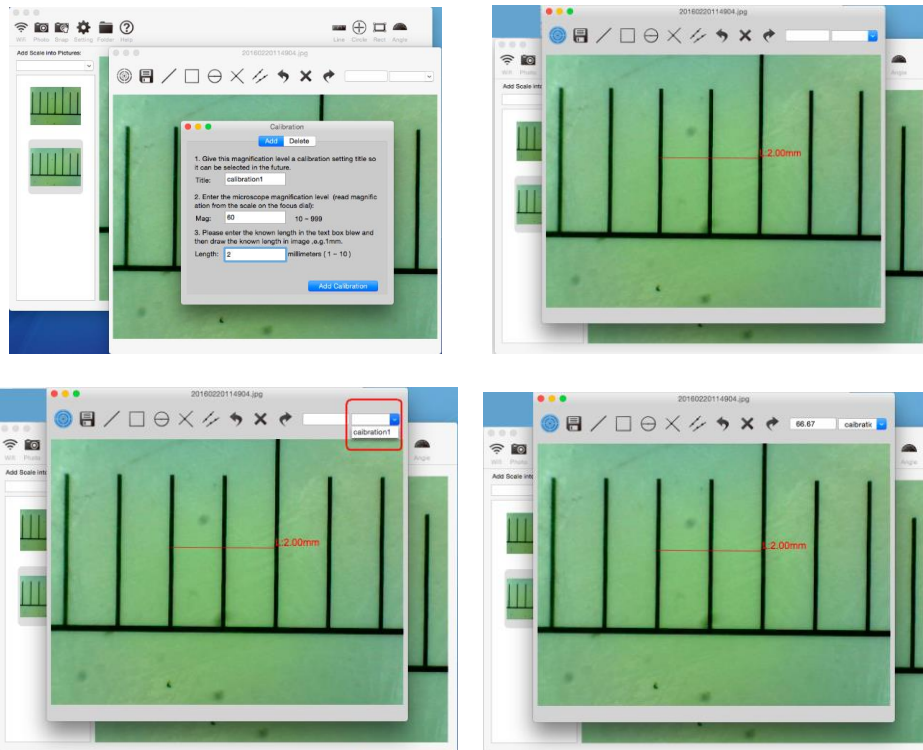Select Magnification on the left-top of window. Click the  icon button; you will see one line on the video preview window. If you select the left of the line and drag the mouse, you can move the line. Select the right of the line and drag the mouse, you can increase or decrease the length of the line.

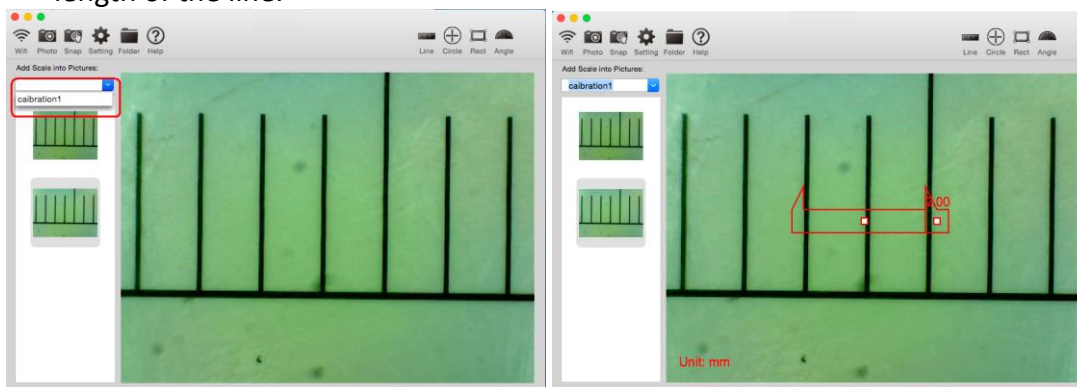
Click  icon button, you will see one circle on the video preview window. If you select center points of the circle and drag the mouse, you can move the circle. Select the right OR top point of the circle and drag the mouse, you can increase or decrease the size of the circle.

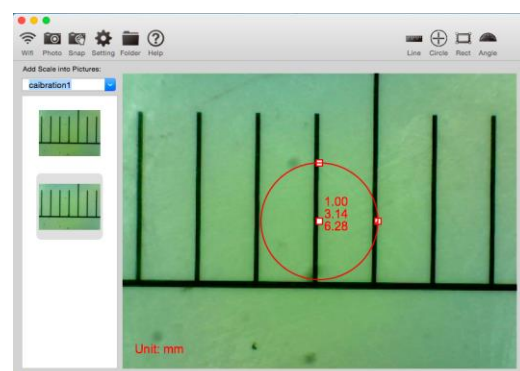
Click **Rect** icon button, the rectangle will display in the video preview window.

Select center point of the rectangle and drag the mouse, you can move the rectangle. Select right-bottom point of the rectangle and drag the mouse, you can change the area of the rectangle.

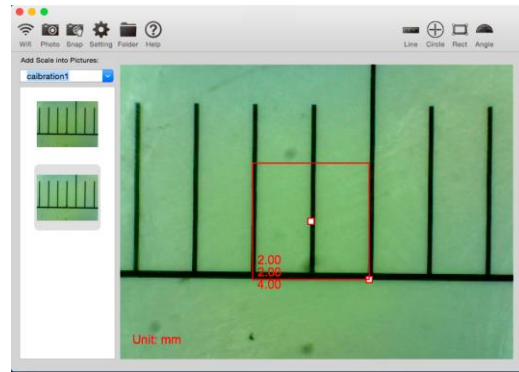
Click **Angle** icon button, the angle will display in the video preview window. Select right OR top point and drag the mouse, you can change the angle;

select corner point of the angle and drag the mouse, you can move the angle.

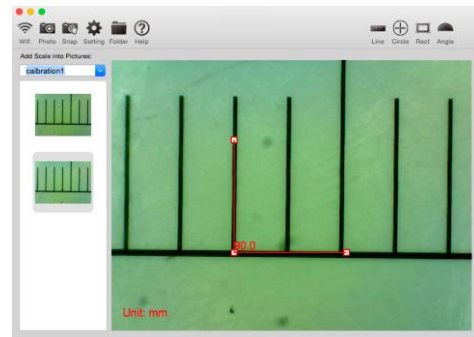Below, the list of the keys available after double click on the image captured

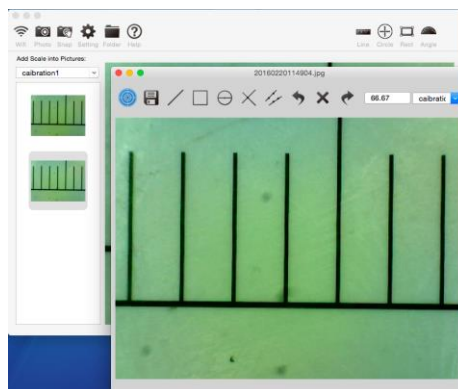Save as: Rename the file and save in different path.

Line: allows you to view the linear distance between two points. Simply left click and drag to the desired length, and click again to finish.

Rectangle: The radius of circle allows you to view the width, height, and area of a rectangle.

Diameter Circle: The diameter of a circle allows you to view the circumference, area, and radius of a circle. Simply click and extend out to the desired diameter.

Three Points Angle: The three point angle allows you to view of an angle.

Four Points Angle: The four point angle allows you to view of an angle by selecting four points.

Undo: Undo the last action.

Delete: Delete.

Redo: Redo the previous action.

13.10 Working at Wired USBMode

2. Install the “**USBMic-Fi**” (for **Windows OS**) or “**USBMic-Fi**” (for **Mac OS**) software; the download is available from www.mic-fi.com
3. Connect to PC’s USB port with USB cable; the power LED  indicator will light, as well as charge indicator LED  (Battery is being charged)
4. Launch “USBMic-Fi”. The illumination LED lights up as soon as the application software is running.


13.11 “USBMic-Fi” on Windows

- Windows XP SP2 or later
- Available USB 2.0 or later port


13.11.1 Running “USBMic-Fi”

Note: Before start USBMic-Fi program, please make sure USB digital microscope has been connected to PC via USB line and Power on.


After installation of USBMic-Fi.exe, double-Click USBMic-Fi icon  on desktop , the USBMic-Fi main screen will be shown:

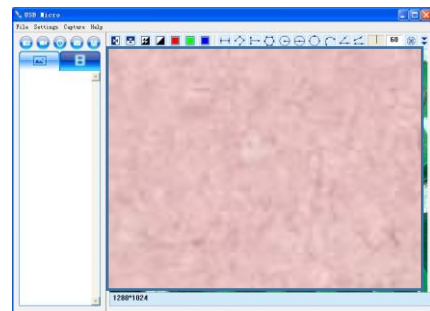If USB digital microscope has not been found while opening the USBMic-Fi program, the message “please connect USB Digital Microscope” will be shown:

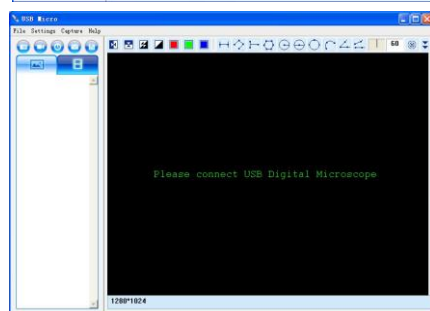13.11.2 Menu

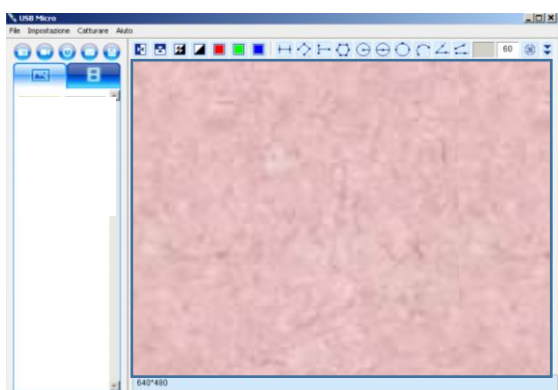Take a Picture


Take/Stop Video

Take time-lapsed pictures


Email the selected photos

Save the selected photo or video file



or

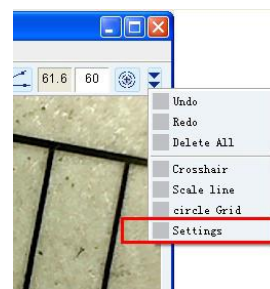Select a photo or video in Photo or Video folder; click “File” on menu bar or right key of mouse to pup-up the drop-down menu for operating it.

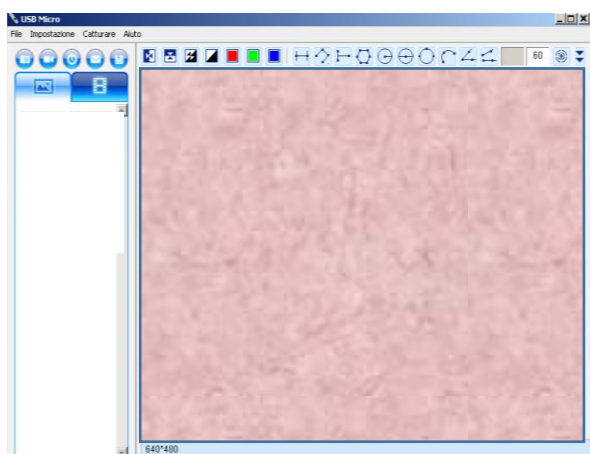- **Open:** Opens the selected picture or video to be viewed
- **Save as**  : Save your files into the desired folder and file format
- **Save all:** Save all your files in the photo or video folder into the desired folder
- **Delete:** delete the selected image or video. You will be prompted to confirm your decision.
- **Delete All:** delete all the images or videos in the photo or video folder. You will be prompted to confirm your decision.
- **Print:** Print the selected photo
- **Email**  : Send the selected images through your default email account.

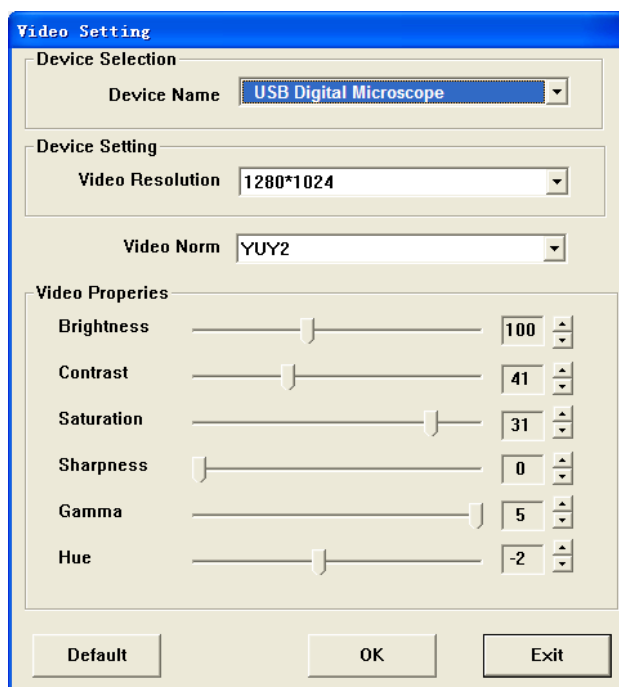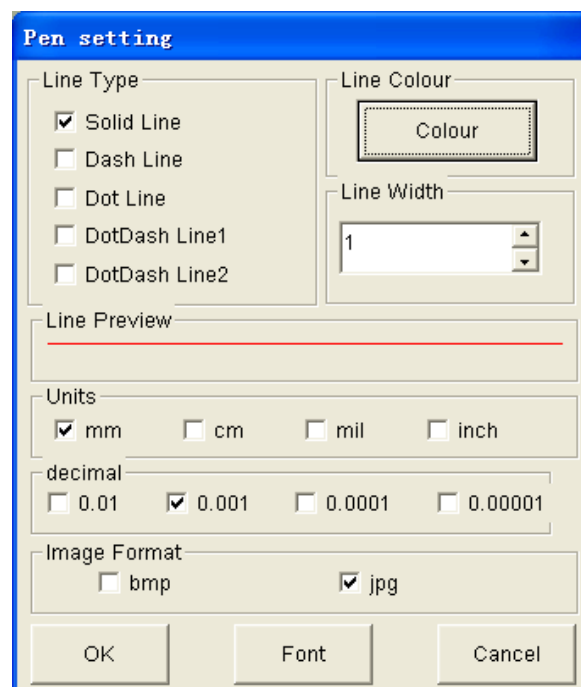Video Setting**Pen Setting**

Crosshair / Scale Line / Circle Grid: Add in Crosshair / Scale Line / Circle Grid on preview area: click “Crosshair / Scale Line / Circle Grid” on the drop-down menu as below.


Snapshot Enable and Snapshot Disable: the default is enable snapshot.

Suggestion: If you need more performance of the video preview, please disable it, or maybe there is video delay.

Change Language

The language selection allows you to select the USBMic-Fi interface language.

Take picture 📷: You can take image by clicking the snapshot Button or press the button 📷 on the USB DigitalMicroscope

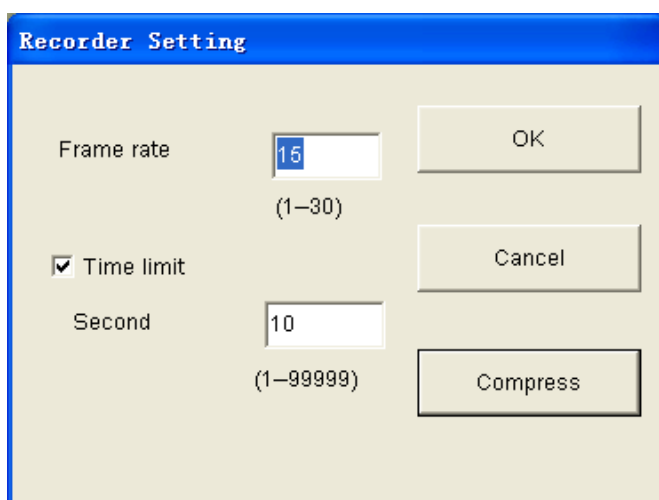**Take Video/Stop Video:**

1. **Frame Rate:** it is the number of frames or images that are projected or displayed per second. The number range is 1~30
2. **Timer Limit:** it is the time of video recording to avoid video file size exceed to the space of disk. The number range is 1~99999
3. **Compress:** click it and you can select the video compressed format to reduce video size. It loads automatically compress programs according to your PC system. The video quality may be varied with different compress program.

Capture Setting

Start Time Year Month Day
 Hour Minute second

Interval Minute second

Capture Frame

File Size MB

Completed Frame

Take time lapsed images:

Press the time-lapsed button  on the active preview window

About USBMic-Fi: you can check your present version Video Process Tool Bar

-  **HorizontalMirror**
-  **VerticalMirror**
-  **Negative**
-  **Binary Segmentation**
-  **Monochrome display (Red)**
-  **Monochrome display (Green)**
-  **Monochrome display (Blue)**
-  **Calibration**

13.11.4 Calibration on Windows

See paragraph 13.8.5

13.12 “USBMic-Fi” on MAC OS

- MAC OS X 10.8 or later
- QuickTime 7.0 or later

Please drag “Mic-Fi” Software to “Applications” Folder to complete installation.

13.12.1 Running “USBMic-Fi”


Install the “USBMic-Fi” software from www.mic-fi.com.


Once the USBMic-Fi installation is completed, please plug the microscope into an available USB port on your Mac OSX PC, and turn the microscope on with the switch on the side. The system will automatically search and install the driver. The software

installation is completed, and you may see a transmission from the microscope on your screen as below.

If not, please read the following instructions. After installation of USBMic-Fi.dmg, double-Click USBMic-Fi  icon , the USBMic-Fi main screen will be shown:

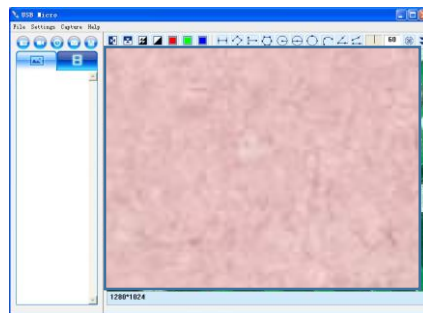If images from the microscope camera are not being transmitted and display is still showing the “Please connect USB Microscope, then click preview...” on screen, please try to change low resolution “Controls →Format-> 320X240” like below.

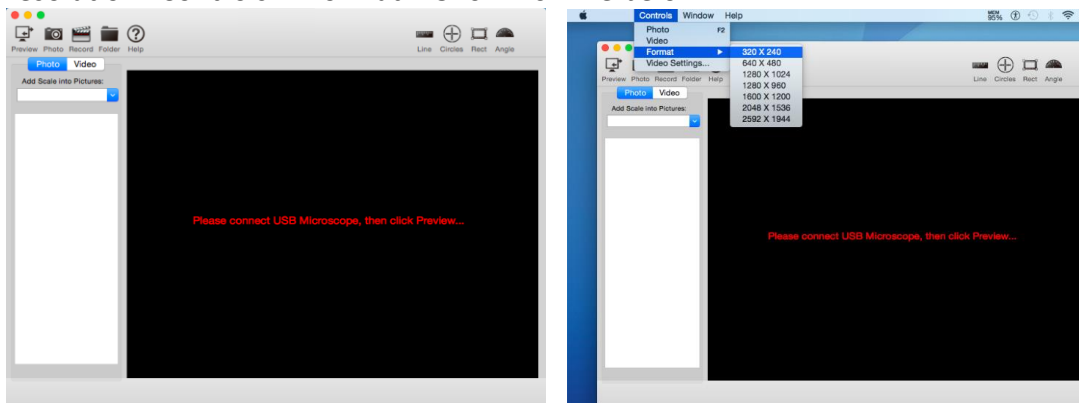13.12.2 *Setting*

Note: When you are planning to use the microscope in the future, be sure to plug it in and turn it on before opening Mic-Fi

Video Resolution Setting

Click “Controls->Format” to select your desired resolution as below screen. The default resolution is 640X480. If your Mic-Fi screen freezes on an image, it may be because capture resolution has been set too high. Set it back to Medium and the video display will return.

Click the  icon button on the Tool bar of the main screen to capture an image. The Captured image will be save in the folder (../users/local user name/Mic-Fi/Image). The filename of the captured image is a unique time stamp with .jpg file. The format of the time stamp filename is: YYYYMMDDhhmmss.jpg : YYYY is the year, MM is month, DD is date, hh is hour, mm is minutes, and ss is second.

Each captured image will be added to the image scroll list. Double-click a selected image thumbnail in the image scroll list to open the image for editing.

Delete a Captured image: Select the captured image in the image scroll list -> press “Del” Key of keyboard to delete it.

Click the  icon button to start recording video and stop recording video on the Tool bar of the main screen as below. The Video filename of the time stamp is YYYYMMDDhhmmss.mov: YYYY is Year, MM is Month, DD is Date, hh is Hour, mm is Minutes, and ss is Second.

The Created Video is saved in the folder (../users/local user name/Mic-Fi/Video). Double-click a selected Video for display or press “Del” key to delete it.

Click the  icon button to View the captured images and video path.

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