**MOPEDIA®**

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---

**MANUALE D'ISTRUZIONE**

## INDICE

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 **MORETTI S.P.A.**

Via Bruxelles, 3 - Meleto 52022 Cavriglia (Arezzo) Tel. +39 055 96 21 11


**www.morettispa.com** email: **info@morettispa.com****MADE IN P.R.C.**

**MORETTI S.P.A.**

Via Bruxelles, 3 - Meleto  
52022 Cavriglia (Arezzo)

Tel. +39 055 96 21 11  
Fax. +39 055 96 21 200

[www.morettispa.com](http://www.morettispa.com)  
[info@morettispa.com](mailto:info@morettispa.com)

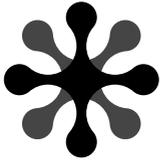**MOPEDIA®**

MA RP922-RP924-RP926 01A ENG\_01+2019

**PEDAL EXERCISER**

---

**INSTRUCTION MANUAL**

## INDEX

<b>1. CODES</b> .....	pag.3
<b>2. INTRODUCTION</b> .....	pag.3
<b>3. INTENDED USE</b> .....	pag.3
<b>4. DECLARATION OF EC CONFORMITY</b> .....	pag.3
<b>5. WARNINGS</b> .....	pag.4
5.1 General warnings .....	pag.4
5.2 Warnings for use.....	pag.4
<b>6. SYMBOLS</b> .....	pag.4
<b>7. GENERAL DESCRIPTION</b> .....	pag.4
7.1 Description of parts RP922 - CROSSY FIX .....	pag.5
7.2 Description of parts RP924 - CROSSY HOLDY.....	pag.5
7.3 Description of parts RP926 - CROSSY TOP .....	pag.5
<b>8. BEFORE USE</b> .....	pag.6
<b>9. ASSEMBLY</b> .....	pag.6
<b>10. INSTRUCTION FOR USE</b> .....	pag.7
10.1 Use with legs .....	pag.8
10.2 Use with arms .....	pag.8
<b>11. MAINTENANCE</b> .....	pag.8
<b>12. CLEANING AND DISINFECTION</b> .....	pag.8
12.1 Cleaning.....	pag.8
12.2 Disinfection .....	pag.9
<b>13. CONDITIONS OF DISPOSAL</b> .....	pag.9
<b>14. SPARE PARTS AND ACCESSORIES</b> .....	pag.9
<b>15. TECHNICAL FEATURES</b> .....	pag.9
15.1 Dimensions and weight RP922 - CROSSY FIX.....	pag.9
15.2 Dimensions and weight RP924 - CROSSY HOLDY .....	pag.10
15.3 Dimensions and weight RP926 - CROSSY TOP.....	pag.10
<b>16. WARRANTY</b> .....	pag.10
<b>17. REPAIRING</b> .....	pag.10
<b>18. SPARE PARTS</b> .....	pag.11
<b>19. EXEMPT CLAUSES</b> .....	pag.11

**CE** Class I medical device

according Directive 93/42/EEC and subsequent amendments

## 1. CODES

- RP922**    CROSSY FIX - Pedal exerciser in painted steel  
**RP924**    CROSSY FOLDY - Folding pedal exerciser in painted steel  
**RP926**    CROSSY TOP - Folding pedal exerciser in painted steel, with display

## 2. INTRODUCTION

Thank you for choosing a **MOPEDIA** pedal exerciser by Moretti S.p.A. The Moretti pedal exercisers are designed and built to meet all your demands for a practical, correct and safe use. This user manual provides some suggestions as to how correctly use the device you have chosen and offers a lots of valuable advice for your safety. Please read carefully through the manual before using the pedal exerciser. Should you have any queries, please contact your retailer for advice and assistance.

**NOTE** Check that all parts of the product have not been damaged during the shipment. In case of damage do not use the product and contact the retailer for further instructions.

### WARNING!

- Do not use the product for a purpose not indicated in this manual.
- Moretti S.P.A. declines all responsibilities for any consequences resulting from an incorrect use of this product and from unauthorized alteration to the frame of the product.
- The manufacturer reserves the right to change the information contained in this document without previous notice.

## 3. INTENDED USE

The pedal exerciser is a medical device indicated for people with mobility difficulties and it is indicated to reestablish coordination and the muscle mass of arms and legs.

## 4. DECLARATION OF EC CONFORMITY

MORETTI S.p.A. hereby declares, under its own exclusive responsibility, that the products manufactured and sold by MORETTI S.p.A. and belonging to the PEDAL EXERCISERS category are compliant with the applicable provisions outlined in the Medical Devices Directive 93/42/EEC. To this aim MORETTI S.p.A., under its own responsibility, guarantees and declares what follows:

1. The devices in question satisfy the essential requirements of Annex I of directive 93/42/EEC as prescribed by Annex VII of the mentioned directive. The complete list of the devices in question is available in annex A to this declaration;
2. The devices in question **ARE NOT MEASUREMENT TOOLS**;
3. The devices in question **ARE NOT INTENDED FOR CLINICAL INVESTIGATION**;
4. The devices in question are sold in **NON STERILE PACKING**;
5. The devices in question should be considered as belonging to class I;
6. **DO NOT USE/INSTALL** the devices in question for any other purpose than that one expressly indicated by MORETTI S.P.A.;
7. Moretti S.P.A. maintains and submits the technical documentation certifying the conformity of the products in question with the directive 93/42/EEC to the proper authorities.

## 5. WARNINGS

### 5.1 General warnings

#### WARNING!

Always pay close attention to the presence of moving parts that could cause limb entrapment and personal injuries.

- For a proper use of the device always consult your doctor or therapist.
- Keep the packed product away from any heat source as the packaging is made of carton.
- The useful life of the accessory is determined by wear of parts that cannot be repaired and/or replaced.
- Be careful when children are nearby.

### 5.2 Warnings for use

- All the legs must lean against the surface in every moment during the use, in order to guarantee the correct balance of the pedal exerciser and the maximum safety.
- Do not straighten and reuse a deformed pedal exerciser. Risk of accident!
- Unauthorized changes to the pedal exercisers or use of non-original spare parts may change the product frame and create conditions of danger, which may cause serious damages and invalidate the warranty.
- The product is subject to particular efforts that depend on the quality of use and from the weight of the person. For this reason we suggest to check the integrity of the product every 3/4 months, especially if it is subjected to impacts or unnatural stress for its use.

## 6. SYMBOLS

Product code

CE mark

Producer

Month and year of production

Read the instruction manual

Conditions of disposal


Disposal

## 7. GENERAL DESCRIPTION

Every CROSSY pedal exerciser of the line MopediA by Moretti is equipped with four non-slip tips, pedals with adjustable foot straps, an handle for the intensity regulation and display (only for the RP926).


### 7.1 Description of parts RP922 - CROSSY FIX

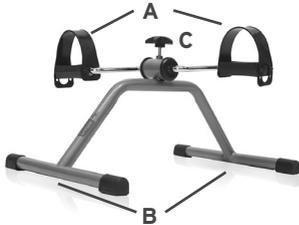List of parts

- A - Pedals with foot straps
- B - Support legs and non-slip tips
- C - Handle for brake regulation (effort intensity)


### 7.2 Description of parts RP924 - CROSSY FOLDY

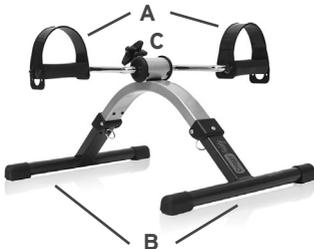List of parts

- A - Pedals with foot straps
- B - Support legs and non-slip tips
- C - Handle for brake regulation (effort intensity)

### 7.3 Description of parts RP926 - CROSSY TOP

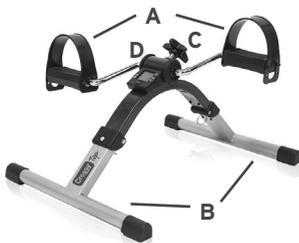List of parts

- A - Pedals with foot straps
- B - Support legs and non-slip tips
- C - Handle for brake regulation (effort intensity)
- D - Display

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## 8. BEFORE USE

- Always verify the state of wear of the mechanical parts, in order to guarantee a safe use of the device, for people and things.
- Verify the correct assembly of the device. Pay particular attention to the blocking of the pedals and to insertion of the tips on the support legs.
- Verify the correct functioning of the brake for the regulation of the intensity effort.

### WARNING!

Always consult your doctor or your therapist before doing every rehabilitative therapy.

## 9. ASSEMBLY

The pedal exerciser is delivered completely or partially assembled and it does not need screws or wrench to be assembled. Before using the pedal exerciser for the first time, verify that there are not damaged parts: if there are, immediately substitute the damaged parts.

Fig.1 - RP922  
Images for illustrative purpose

Fig.2 - RP924, RP926

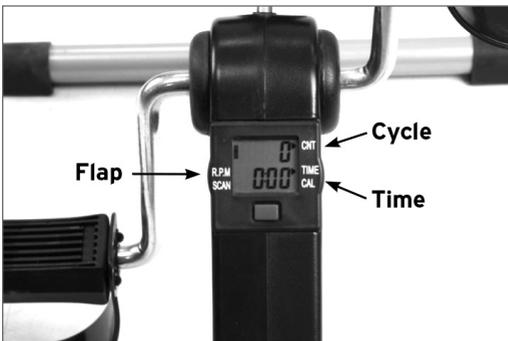
Fig.3  
Images for illustrative purpose

Fig.4

1. Open the box and extract the product;
2. Insert the brake handle in its compartment (Fig.1);
3. Open the pedal exerciser pulling the legs outwards until the click of the snap fastener (only for RP924 and RP926) (Fig.2);
4. To turn on the display the first time (only for RP926) extract the flap positioned to left of the screen, then start pedaling and the display will automatically turn on (Fig.3);
5. The display shows the number of cycles and the running time. To reset it simply hold down the red button below the screen for a few seconds. The display automatically turns off after a short period of inactivity (only for RP926) (Fig.3);
6. To close the pedal exerciser pull outwards the ring of the snap fastener (figure 4) and close the legs inwards (only for RP924 and RP926) (Fig.2).

## 10. INSTRUCTION FOR USE

### WARNING!

The pedal exerciser require regular maintenance and must be daily checked to work in a proper way. We suggest to immediately substitute the worn or broken parts that may cause damages to the patient. It is recommended that an expert shows how to use the pedal exerciser and that the patient repeat the procedure under his/her supervision. For a correct rehabilitation it is necessary to maintain a correct posture. Before beginning the activity, consult your doctor/therapist.

- Ensure that the screws are all completely clamped before use;
- Place the product in a stable way and in a flat surface;
- The positioning of the pedal exerciser in relation to the height of the seat must permit the complete distension of the leg without having to incline the pelvis to the right and to the left;
- Start the activity only when the chest is almost straight;
- To gradually effectuate the rehabilitation, the pedal exerciser is equipped with brake for the regulation of the effort intensity;
- Do not take the feet off the pedals while the pedal exerciser is in fuction;
- At the end of the use, turn clockwise the handle for the brakes regulation, to increase the power and avoid that the pedals freely rotate causing possible injuries.

### WARNING!

The non-compliance of the previous recommendations may cause injuries or serious health problems.

### 10.1 Use with legs

1. Place the pedal exerciser in a flat and non slip surface;
2. Sit down in front of the device at a distance that allows to easily insert the feet on the pedals. **Always place the right foot on the pedal labelled with the letter "R" and the left foot on the pedal labelled with the letter "L".**
3. Adjust the foot straps as necessary;
4. Adjust the brake hardness following the doctor or therapist's recommendations.

**NOTE** To increase the effort intensity turn clockwise the brake handle. To decrease the intensity turn counterclockwise.

5. Activate the device doing a rotatory movement.

### 10.2 Use with arms

1. Place the pedal exerciser in a flat and non slip surface;
2. Sit down in front of the device at a distance that allows to easily insert the feet on the pedals. **Always place the right foot on the pedal labelled with the letter "R" and the left foot on the pedal labelled with the letter "L".**
3. Regulate the hardness of the brake following your doctor or therapist's recommendations.

**NOTE** To increase the effort intensity turn clockwise the brake handle. To decrease the intensity turn counterclockwise.

4. Firmly grab the pedals with the hands and activate the device doing a rotatory movement.

## 11. MAINTENANCE

The MOPEDIA devices by MORETTI are checked carefully and supplied with EC mark, once launched on the market. For patient and doctor's safety we recommend to let check the device by the manufacturer or by laboratory authorized once a year.

For periodic inspection is intended a check of the following parts:

- The supporting framework of the pedal exerciser and its connecting parts;
- Check the assembly of the pedals, of the tips and the functioning of the brake for the regulation of the effort intensity;
- Check that all the bolts are tightened;
- Verify all the welding points;
- Lubrication: the pedal exerciser is thought for a minimal maintenance; however a check every six months and the lubrication should ensure an excellent duration. The pedal exerciser must freely rotate.

In case of repair, use only original spare parts and accessories.

## 12. CLEANING AND DISINFECTION

### 12.1 Cleaning

- Use exclusively a wet cloth and neutral soap. Successively dry well before the use. Do not use direct water jets.
- Use water at a temperature not superior than 30°C.
- Do not use pressure washer or steam.

## 12.2 Disinfection

- If necessary effectuate a disinfection of the product with a common disinfectant cleaner.
- A total disinfection of the product has to be effectuated regularly or when the product it is used by another user. The colour loss of some materials may be caused by particular kind of perspiration, by UV radiations, wear, heat, humidity and it is considered a natural process not covered by warranty.

**NOTE** Absolutely avoid the use of acidic, alkaline or solvent products, such as acetone or diluent.

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### **WARNINGS FOR THE PROPER DISPOSAL OF THE PRODUCT, AS REQUIRED BY THE EUROPEAN DIRECTIVE 2002/96/CE:**

At the end of its service life, the product must not be disposed of together with normal municipal waste. It must be delivered to the municipal sorted waste collection centre or to the resellers, in case they provide for this service. Sorted waste disposal reduces the environmental and health consequences of an improper disposal. Furthermore, the materials which make up the product may be recycled, thus resulting in energy and resource savings. In order to highlight the compulsory sorted disposal of the electro-medical equipment, the product itself bears the barred bin symbol.

### **TREATMENT OF DEAD BATTERIES (ONLY FOR RP926 - CROSSY TOP) (European directive 2006/66/CE):**

The symbol on the product indicates that the batteries must not be treated like common household waste. Make sure the batteries are properly disposed; this prevents the possible environmental and health consequences of an improper disposal. Recycling helps to preserve natural resources. Deliver dead batteries to a proper recycling centre. For more detailed information regarding the disposal of dead batteries or of the product itself, please contact the municipality, the local waste disposal centre or the reseller from who you purchased the device from.

## 14. SPARE PARTS AND ACCESSORIES

For spare parts and accessories, please always refer to the MORETTI master catalogue.

## 15. TECHNICAL FEATURES

Materials: powder coated tubular steel frame.

### 15.1 Dimensions and weight RP922 - CROSSY FIX

External dimensions	54x41x35,5 cm
Minimum height (pedals in horizontal position)	28 cm
Maximum height (pedals in vertical position)	35,5 cm
Weight of the product	2,2 kg

### 15.2 Dimensions and weight RP924 - CROSSY FOLDY

External dimensions	48x37x34 cm
Minimum height (pedals in horizontal position)	27 cm
Maximum height (pedals in vertical position)	34 cm
Weight of the product	2,3 kg

### 15.3 Dimensioni e peso RP926 - CROSSY TOP

External dimensions	49x37x34 cm
Minimum height (pedals in horizontal position)	27 cm
Maximum height (pedals in vertical position)	34 cm
Weight of the product	2,3 kg

## 16. WARRANTY

Every Moretti product is guaranteed for 2 (two) YEARS from the purchase date, except the restrictions specified below. The warranty shall not apply in case of damage caused by improper use, abuse, alterations of the product or any other use not in compliance with the user manual. The proper use of the product is indicated in the user manual.

Moretti S.p.A. shall not be liable for any damages or personal injuries due to incorrect assembly/improper use of the device, not compliant with the user manual. Moretti does not guarantee its products against damages or faults in the following cases: natural disasters, non authorized maintenance or reparations, damages caused by electrical supply problems (if provided), use of parts or components not provided by Moretti, failure in accomplishing the guidelines and instructions for use, non authorized alterations, damages due to the transport (different from Moretti original delivery) or failure in maintaining the device as indicated in the manual. The components subject to wear are not covered by this warranty if the damage is due to the normal use of the product.

## 17. REPAIRING

### Warranty repair

If a Moretti product shows material or manufacturing faults during the warranty period, Moretti S.p.A., together with the customer, will consider whether the product fault is covered by warranty. Moretti S.p.A., at its unquestionable discretion, may replace or repair the item under warranty in a Moretti reseller point or at its own qualified offices. Labour costs for repairing the product may be at Moretti's expense if it is determined the reparation is under warranty. A reparation or replacement does not renew nor extend the warranty.

### Repair not covered by warranty

A product, not covered by warranty, may be returned for reparation only after receiving prior authorization from Moretti customer service. Labour and delivery costs for an out-of-warranty reparation will be fully at the customer's or reseller's expense. Reparations of out-of-warranty products are guaranteed for 6 (six) months, from the day the customer receives the repaired product.

### Non defective devices

The customer is informed if, after examining and testing a returning product, Moretti decides the product is not defective. The product will be sent back to the customer and all the return delivery costs will be at his expense.

## 18. SPARE PARTS

The original Moretti spare parts are guaranteed for 6 (six) months from the day you receive them.

## 19. EXEMPT CLAUSES

Except what clearly specified in this warranty and in accordance with the law, Moretti S.p.A. does not offer its customers any other declaration, warranty or condition, expressed or implicit, including any declaration, warranty or condition of merchantability, of non violation or non interference, of suitability for a specific aim. Moretti S.p.A. does not guarantee that the use of Moretti products will be continuous and faultless. The duration of possible implicit warranties that may be applied in accordance with the law is limited to the warranty period.

Some states or countries do not allow limitations to the duration of an implicit warranty, or the exoneration, or the limitations of accidental or indirect damages related to products for consumers. In such states or countries some exoneration or limitation of this warranty may not be applied to the user. This warranty is subject to modifications without any prior notice.

## WARRANTY CERTIFICATE

Product \_\_\_\_\_

Purchased on (date) \_\_\_\_\_

Retailer \_\_\_\_\_

Address \_\_\_\_\_ Town/city \_\_\_\_\_

Sold to \_\_\_\_\_

Address \_\_\_\_\_ Town/city \_\_\_\_\_

 MORETTI S.P.A.

Via Bruxelles, 3 - Meleto 52022 Cavriglia (Arezzo) Tel. +39 055 96 21 11

[www.morettispa.com](http://www.morettispa.com) email: [info@morettispa.com](mailto:info@morettispa.com)

MADE IN P.R.C.

**MORETTI S.P.A.**

Via Bruxelles, 3 - Meleto  
52022 Cavriglia (Arezzo)

Tel. +39 055 96 21 11  
Fax. +39 055 96 21 200

[www.morettispa.com](http://www.morettispa.com)  
[info@morettispa.com](mailto:info@morettispa.com)